

Stranded Options for LDR Implants



PRE-PLANNED & PRE-LOADED



Innovative stranded options backed by Theragenics' trusted experience and quality.

TheraStrand®

- Braided 90/10 (glycolide/L-lactide) suture containing seeds & spacers
 - Sterile
 - Bio-absorption occurs in 56-70 days
- Variable seed and spacer configurations customized per a physician's plan
 - 2.75 mm, 5.5 mm & 1.0 cm spacer lengths
 - Maximum length of 6.5 cm
- Available in pre-loaded needle trays

TheraSleeve®

- Bioabsorbable 20/80 (glycolide/L-lactide) sleeve containing seeds & spacers
 - Sterile
 - Bio-absorption occurs in 140-180 days
- Variable seed and spacer configurations customized per a physician's plan
 - 2.75 mm, 5.5 mm & 1.0 cm spacer lengths
 - Maximum length of 7.0 cm
- Available in pre-loaded needle trays

Content intended for the USA. Indications may vary in different countries; please either consult the appropriate IFU or your local distributor.



theragenics.com



customerservice@theragenics.com




877.444.7333

Stranded Options for LDR Implants



REAL-TIME STRANDED OPTIONS

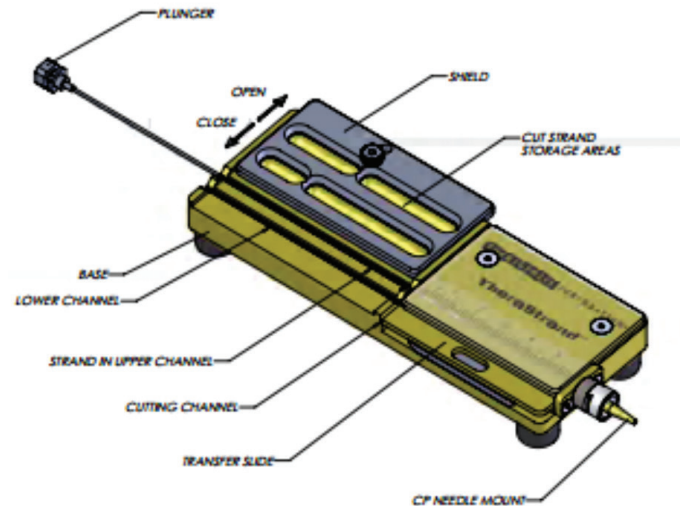

TheraStrand® RT

- Braided 90/10 (glycolide/L-lactide) suture containing seeds & spacers
 - Sterile
 - Bio-absorption occurs in 56-70 days
- Accommodates real-time planning and intra-operative needle loading
- Fixed seed spacing
 - 1 cm spacing center-to-center of seed
 - Maximum strand length of 6 seeds and 6 spacers
 - 2, 3, 4, & 5 fixed seed and spacer strand lengths available upon request
- 2 packaging options
 - 20 strand capacity shielded container
 - 4 strand capacity shielded container suitable for sterile strand well-chamber assay

TheraStrand® RT Cutting Block

- Facilitates sterile TheraStrand RT custom cutting and intra-operative needle loading
 - Reusable
 - High contrast anodized aluminum enhances seed and spacer visibility
 - Precision cutting
- Shielded to reduce handling exposure in accordance with ALARA practices

TheraSeed and AgX100 are the foundation for all LDR brachytherapy product presentations. A manufacturer's Certificate of Calibration accompanies each order; an additional independent seed assay is available upon request.

Content intended for the USA. Indications may vary in different countries; please either consult the appropriate IFU or your local distributor. TheraGenics Corporation, TheraSeed, AgX100 and AnchorSeed are registered trademarks of TheraGenics Corporation. ©2021 TheraGenics Corporation

**For specific product information, please consult your TheraGenics Brachytherapy Specialist,
Customer Service, or our website for the eIFU.**



theragenics.com



customerservice@theragenics.com



877.444.7333