

Brachytherapy solutions to better serve physicians and their patients


TheraSeed[®] & AgX100[®]

Stranded options to meet your LDR brachytherapy implant needs...

Pre-planned & Pre-loaded


TheraStrand[®]

- Braided 90/10 (glycolide/L-lactide) suture containing seeds & spacers
 - Sterile
 - Bio-absorption occurs in 56-70 days
- Variable seed and spacer configurations customized per a physician's plan
 - 2.75 mm, 5.5 mm & 1.0 cm spacer lengths
 - Maximum length of 6.5 cm
- Available in pre-loaded needle trays

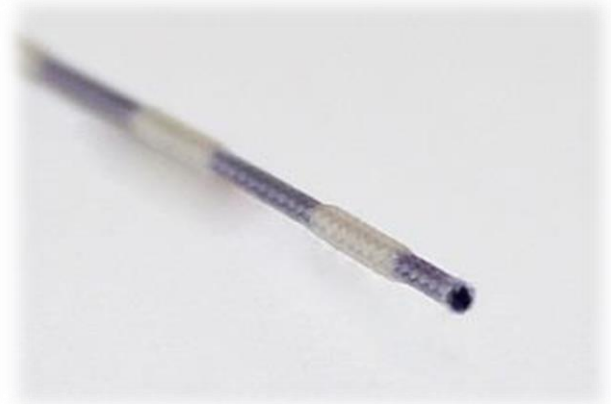

TheraSleeve[®]

- Bioabsorbable 20/80 (glycolide/L-lactide) sleeve containing seeds & spacers
 - Sterile
 - Bio-absorption occurs in 140-180 days
- Variable seed and spacer configurations customized per a physician's plan
 - 2.75 mm, 5.5 mm & 1.0 cm spacer lengths
 - Maximum length of 7.0 cm
- Available in pre-loaded needle trays

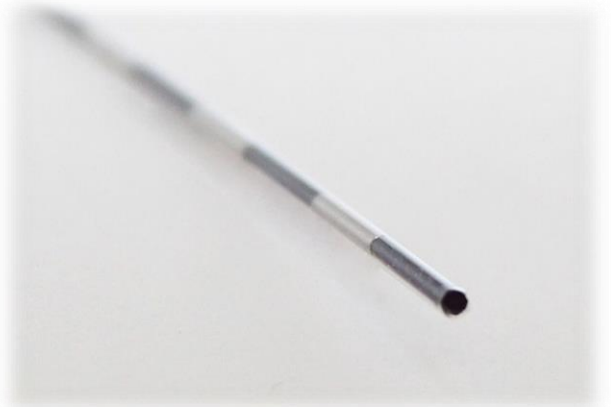

Real-time Stranded Option


TheraStrand® RT

- Braided 90/10 (glycolide/L-lactide) suture containing seeds & spacers
 - Sterile
 - Bio-absorption occurs in 56-70 days
- Accommodates real-time planning and intra-operative needle loading
- Fixed seed spacing
 - 1 cm spacing center to center of seed
 - Maximum strand length of 6 seeds and 6 spacers
 - ❖ 2, 3, 4, & 5 fixed seed and spacer strand lengths available upon request
- 2 packaging options
 - 20 strand capacity shielded container
 - 4 strand capacity shielded container suitable for sterile strand well-chamber assay



TheraStrand® RT Cutting Block

- Facilitates, sterile TheraStrand RT custom cutting and intra-operative needle loading
 - Reusable
 - High contrast anodized stainless steel enhances seed and spacer visibility
 - Precision cutting
- Shielded to reduce handling exposure in accordance with ALARA practices

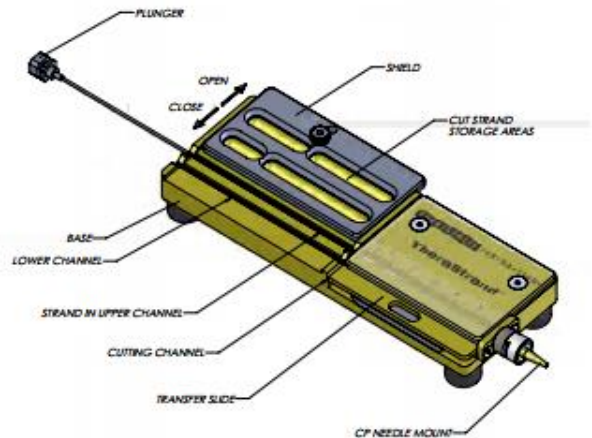

TheraSeed and AgX100 are the foundation for all LDR brachytherapy product presentations. A manufacturers Certificate of Calibration accompanies each order; an additional independent seed assay of 10% is available upon request.

For specific product information and Instructions for Use, please consult your Theragenics Brachytherapy Specialist or Customer Service.

Customer Service

Phone: 877-444-7333

Fax: 800-458-4303

Email: customerservice@theragenics.com

Hours of Operation: Mon - Fri, 8am - 6pm ET

Theragenics Corporation

5203 Bristol Industrial Way

Buford, Georgia 30518

USA