


Catalog No.	Sheath Size	Dilator Length	Wire	Needle Size	Sheath ID	Sheath OD	Tip ID
CLI-230-45	4F 4F x 7cm	11.8cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.066"	0.089"	0.061"
CLI-230-55	5F 5F x 7cm	11.8cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.079"	0.102"	0.074"
CLI-230-06	6F 6F x 7cm	11.8cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.086"	0.109"	0.086"
CLI-230-07	7F 7F x 7cm	12.0cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.099"	0.122"	0.094"
CLI-228-45	4F 4F x 7cm	11.8cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.066"	0.089"	0.061"
CLI-228-55	5F 5F x 7cm	11.8cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.079"	0.102"	0.074"
CLI-228-06	6F 6F x 7cm	11.8cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.086"	0.109"	0.086"
CLI-228-07	7F 7F x 7cm	12.0cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.099"	0.122"	0.094"
CLI-231-45	4F 4F x 11cm	15.8cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.066"	0.089"	0.061"
CLI-231-55	5F 5F x 11cm	15.8cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.079"	0.102"	0.074"
CLI-231-06	6F 6F x 11cm	15.8cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.086"	0.109"	0.086"
CLI-231-07	7F 7F x 11cm	16.0cm	.018" x 45cm NT/Au	21G 21G x 2.5cm	0.099"	0.122"	0.094"
CLI-229-45	4F 4F x 11cm	15.8cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.066"	0.089"	0.061"
CLI-229-55	5F 5F x 11cm	15.8cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.079"	0.102"	0.074"
CLI-229-06	6F 6F x 11cm	15.8cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.086"	0.109"	0.086"
CLI-229-07	7F 7F x 11cm	16.0cm	.018" x 45cm SS/SS	21G 21G x 2.5cm	0.099"	0.122"	0.094"

Quantity: 5 per box

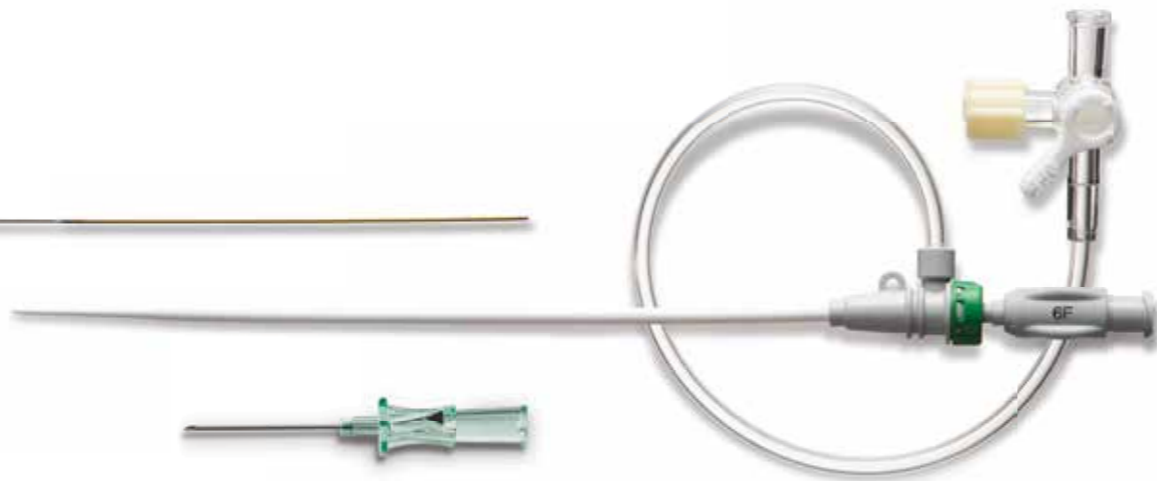

## Precision and Performance

The ELITE HV<sup>®</sup>R transitionless radial introducer sheath is designed for precise access, smooth entry, and consistent performance. Featuring transitionless, hydrophilic sheath properties, combined with a premium guidewire and needle, the ELITE HV<sup>®</sup>R kit provides the tools needed for successful, atraumatic radial access.

**23% Less Insertion Force**  
than a comparable sheath\*

**Hydrophilic Coating**  
for easy insertion and removal

**Transitionless Sheath Tip**  
for smooth entry

\*1:1 simulated data on file

**Theragenics** | [theragenics.com](http://theragenics.com) | [info@theragenicsinterventional.com](mailto:info@theragenicsinterventional.com) | 800.639.2800

Galt Medical Corp. is a subsidiary of Theragenics Corporation

Theragenics Corporation, 5203 Bristol Industrial Way, Buford, Georgia 30518, USA

MKT-200-09 Rev. B


**Theragenics** | [theragenics.com](http://theragenics.com) | [info@theragenicsinterventional.com](mailto:info@theragenicsinterventional.com) | 800.639.2800

Galt Medical Corp. is a subsidiary of Theragenics Corporation

Theragenics Corporation, 5203 Bristol Industrial Way, Buford, Georgia 30518, USA

MKT-200-09 Rev. B

Name: 21THE020\_Elite\_HV\_R\_SS\_MKT\_200\_09\_Rev\_A\_17x11\_R4\_M  
Date: 04/18/22 Version/Misc. N/A  
Notes: N/A

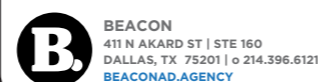


AE	CD	CW	PMS
AD	ST	PR	PMS

Trim: 17"W x 11"H  
Live: 17"W x 11"H  
Folded: 8.5"W x 11"H  
Bleed: .125" (to all sides)  
Scale: 100% Res: 300dpi  
Mech Size: 17"W x 11"H

TRIM

## Radial Artery Access Sheath

Sheath/dilator assembly accommodates an .018" guidewire and is available in sizes from 4F to 7F and in 7cm and 11cm lengths.

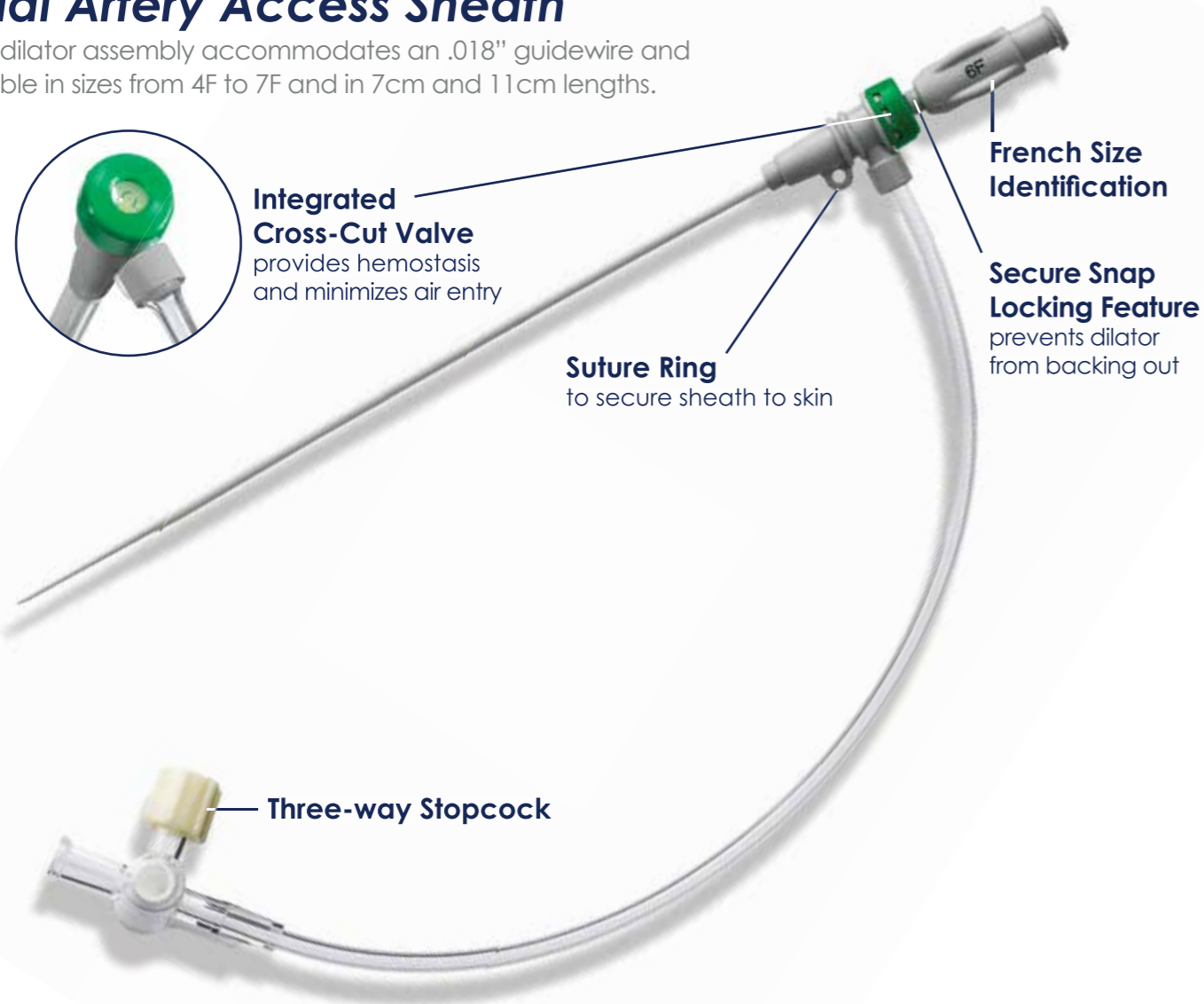

**Integrated Cross-Cut Valve** provides hemostasis and minimizes air entry


**French Size Identification**

**Secure Snap Locking Feature** prevents dilator from backing out

**Suture Ring** to secure sheath to skin

**Three-way Stopcock**

**Sheath Tip ID/Sheath OD Specifications**

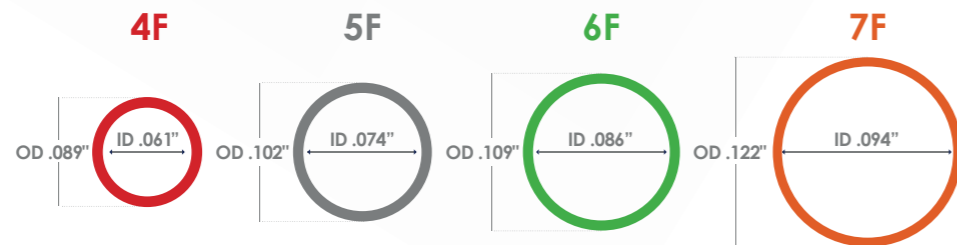

**Theragenics** | [theragenics.com](http://theragenics.com) | [info@theragenicsinterventional.com](mailto:info@theragenicsinterventional.com) | 800.639.2800


Galt Medical Corp. is a subsidiary of Theragenics Corporation

Theragenics Corporation, 5203 Bristol Industrial Way, Buford, Georgia 30518, USA

MKT-200-09 Rev. B

## Guidewire

Kits include an .018" x 45cm guidewire with nitinol core and gold radiopaque tip (shown below) or stainless steel tip and core. Separate sterile guidewires in additional lengths are also available.

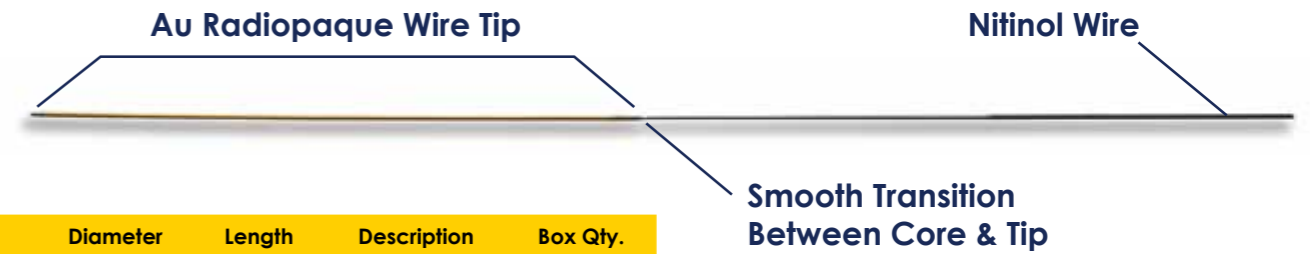

Catalog No.	Diameter	Length	Description	Box Qty.
SGW-056-01	.018"	40cm	NT/Au	5
SGW-056-04	.018"	45cm	NT/Au	5
SGW-015-02	.018"	40cm	SS/SS	10
SGW-015-08	.018"	45cm	SS/SS	10

## 21G Needle

Kits include a premium 21G x 2.5cm needle, designed for precision and featuring a bevel indicator on the hub. Separate sterile needles, including X Sharp®, are also available.



Catalog No.	Needle Size	Description	Box Qty.
NDL-100-32	21G x 2.5cm	Smooth	25
NDL-100-17	21G x 4cm	Smooth	25
NDL-107-14	21G x 4cm	Smooth, X Sharp	25

**Theragenics** | [theragenics.com](http://theragenics.com) | [info@theragenicsinterventional.com](mailto:info@theragenicsinterventional.com) | 800.639.2800

Galt Medical Corp. is a subsidiary of Theragenics Corporation

Theragenics Corporation, 5203 Bristol Industrial Way, Buford, Georgia 30518, USA

MKT-200-09 Rev. B

Name: 21THE020\_Elite\_HV\_R\_SS\_MKT\_200\_09\_Rev\_A 17x11\_R4\_M  
Date: 04/18/22 Version/Misc. N/A  
Notes: N/A

**B.** BEACON  
411 N AKARD ST | STE 160  
DALLAS, TX 75201 | o 214.396.6121  
BEACONAD.AGENCY

AE	CD	CW	PMS
AD	ST	PR	PMS

Trim: 17"W x 11"H  
Live: 17"W x 11"H  
Folded: 8.5"W x 11"H  
Bleed: .125" (to all sides)  
Scale: 100% Res: 300dpi  
Mech Size: 17"W x 11"H